INNOVATIVE DESIGNS  
for independent living™


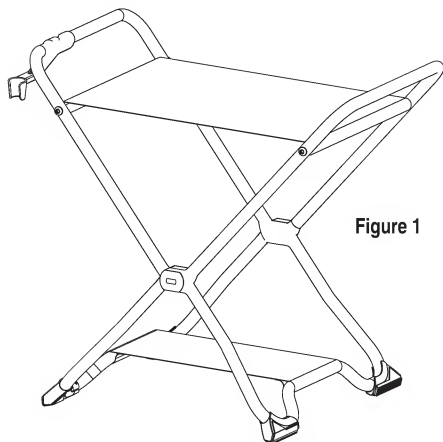
Figure 1


# DN7026 HOME CARE® BY MOEN® FOLDING MESH SHOWER SEAT INSTRUCTIONS

INSTRUCCIONES PARA EL ASIENTO DE MALLA PLEGABLE PARA REGADERA HOME CARE DE MOEN  
HOME CARE DE MOEN DIRECTIVES D'INSTALLATION POUR SIÈGE DE DOUCHE PLIANT EN FILET

To Remove the Handheld Shower Holder/Para quitar el soporte para regadera de mano/Pour enlever le porte-douche à main

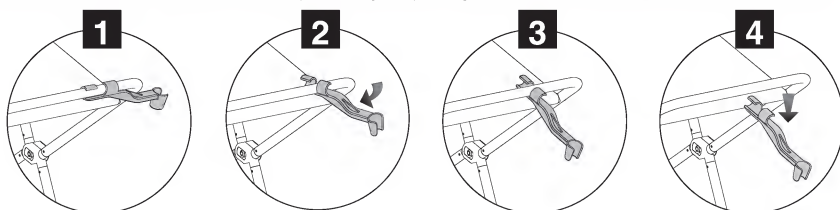To Attach the Handheld Shower Holder/Para fijar el soporte para regadera de mano/Pour fixer le porte-douche à main

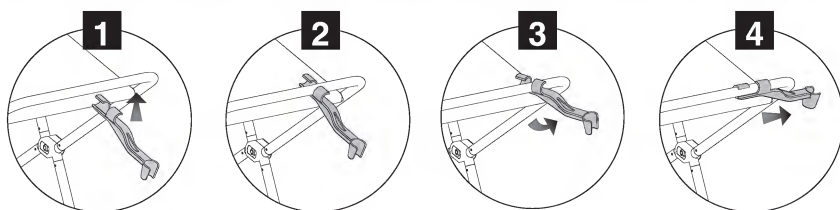Thank you for purchasing this Home Care® by Moen® product. Home Care by Moen is a leader in offering top-quality products that are attractive and innovative, with performance features like a 250 lb. weight capacity that exceeds ADA requirements.

The entire Home Care by Moen family of products offers people new options for independent living by providing outstanding comfort, convenience and style.

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## Instructions

### NO TOOLS REQUIRED!

- 1). Remove Folding Mesh Shower Seat from plastic bag.
- 2). Grabbing the arms of the seat firmly, pull them in opposite directions until the seat is fully opened.
- 3). Slide the Handheld Shower Holder to the end of the armrest and rotate the holder to the upright position (see Figure 1).
- 4). Insert a handheld shower into the holder (not included).

### To Remove the Handheld Shower Holder

- 1). Rotate the Handheld Shower Holder so that it is level with the arm of the Folding Mesh Shower Seat.
- 2). Pull the Handheld Shower Holder toward you until it is in the upright position and the clips are released from the arm of the Folding Mesh Shower Seat.
- 3). Begin to pull the Hand Held Shower Holder away from the arm of the Folding Mesh Shower Seat.
- 4). The Handheld Shower Holder is now released.

### To Attach the Handheld Shower Holder

- 1). Align the Handheld Shower Holder with the arm of the Folding Mesh Shower Seat to which you would like to attach your handheld shower.
- 2). Slide the Handheld Shower Holder clips so that the arm of the Folding Mesh Shower Seat is between the clips.
- 3). Rotate the Handheld Shower Holder toward the arm until the clips snap into place.
- 4). Slide the Handheld Shower Holder forward; be sure it is in the upright position and insert a handheld shower (not included).

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

## Components Included/Componentes incluidos/Pièces incluses

1 Folding Mesh Shower Seat  
1 asiento de malla plegable  
1 siège de douche pliant en filet

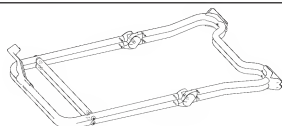1 Handheld Shower Holder  
1 soporte para regadera de mano  
1 porte-douche à main

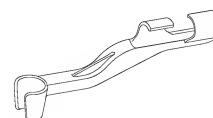INNOVATIVE DESIGNS  
for independent living™

**WARNING:**

TO PREVENT SERIOUS OR FATAL INJURY: FOR USE ONLY BY PERSONS WEIGHING LESS THAN 250 POUNDS. Consult a healthcare professional for advice and instructions specific to your needs before use. Contact CSI immediately at 1-800-882-0116 if any components are missing or damaged. Install in accordance with CSI's instructions. Do not modify this product in any way. FOR QUESTIONS OR ASSISTANCE WITH INSTALLATION CONTACT CSI AT 1-800-882-0116.

**LIMITED LIFETIME WARRANTY**

If this product fails due to a defect in materials or workmanship at any time during the life of the product, Creative Specialties International (CSI), a division of Moen Incorporated, will replace it free of charge, postage-paid. Simply contact CSI at 1-800-882-0116 for replacement information. This warranty does not cover products which have been abused, altered, damaged, misused, cut or worn. TO THE EXTENT PERMITTED BY LAW, CSI DISCLAIMS ALL OTHER IMPLIED OR EXPRESS WARRANTIES INCLUDING ALL WARRANTIES OF MERCHANTABILITY AND/OR FITNESS FOR A PARTICULAR PURPOSE.

**LIMITATION OF LIABILITY:** THIS WARRANTY IS YOUR SOLE REMEDY AND CSI SHALL NOT BE LIABLE FOR ANY DAMAGES, WHETHER DIRECT, INDIRECT, INCIDENTAL, SPECIAL, CONSEQUENTIAL, EXEMPLARY, OR OTHERWISE, INCLUDING LOST REVENUES AND LOST PROFITS, ARISING OUT OF ANY THEORY OR RECOVERY, INCLUDING STATUTORY, CONTRACT OR TORT. NOTWITHSTANDING THE TERM OF ANY LIMITED OR IMPLIED WARRANTY, OR IN THE EVENT THAT ANY LIMITED WARRANTY FAILS OF ITS ESSENTIAL PURPOSE, IN NO EVENT WILL CSI'S ENTIRE LIABILITY EXCEED THE PURCHASE PRICE OF THIS PRODUCT. SOME STATES AND PROVINCES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES SO THE ABOVE LIMITATIONS OR EXCLUSIONS MAY NOT APPLY TO YOU. THIS LIMITED WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS, AND YOU MAY ALSO HAVE OTHER RIGHTS, WHICH VARY FROM STATE TO STATE AND PROVINCE TO PROVINCE.

For additional warranty information, visit [www.homecare.moen.com](http://www.homecare.moen.com)

Home Care by Moen Folding Mesh Shower Seat meets or exceeds the applicable requirements of the following standards/guidelines:

 The Americans with Disabilities Act (ADA)

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Creative Specialties is a division of Moen Incorporated. © 2006
Home Care is a registered trademark of Home Care Tim-BR-MART, licensed to Moen Incorporated.
Home Care and Creative Specialties are a registered trademarks of Moen Incorporated.